Long Term Care products designed to exceed resident requirements for comfort and functionality and surpass customer expectations for high quality, selection and value.

Our broad line includes standard Long Term Care & Bariatric Beds, Resident Room Furniture that is both attractive and practical, Pressure Prevention products that meet a wide range of clinical demands, Multi-Position Recliners, Wheelchairs, Bathing Systems, Bath Safety, Respiratory, and Mobility aids just to name a few.

We strive to be your sole source supplier for all your facilities resident care needs.





## **CONTENTS**

Beds	1
Surfaces	3
Positioning Chairs	12
Patient Room	13
Patient Transfer	17
Bathing	23
Furniture	25
Disposables & Infection Control	35
Fautionant & Diagnastics	27

### **BEDS**

#### Joerns® Healthcare Beds

- Optional roll-at-any-height mobility.
- Come standard with a hand pendant, wall bumper and mattress keepers included.
- Multiple rail, assist and panel options.
- Optional staff control with pendant lock-out.

### Joerns® EasyCare® Healthcare Bed

- Travel range of 7" to 30"
- Standard or AdvanceCare functionality
- Optional integrated width (at time or order) or CareWide width extension available; increases bed width 39" or 42"
- Ergonomic and Intuitive User Controls
- Safe Working Load: 600lbs



### Joerns® UltraCare XT® Healthcare Bed

- Travel range of 7" to 30"
- UltraLock locking casters are available at the foot or head and foot (at time of order)
- QuickCare improves staff efficiency with fast actuators that move patients into position more quickly than other beds
- Ergonomic and Intuitive User Controls
- UltraWide width extension available; increases bed width to 39" or 42"
- Safe Working Load: 500lbs.



### Joerns® Care 100® Healthcare Bed

- Travel range of 8 ¾" to 26"
- CareWide width extension available; increases bed width 39" or 42"
- Ergonomic and Intuitive User Controls
- Safe Working Load: 450lbs.



Permobil is pleased to offer a new bed option to meet a wider variety of patient needs. The new Tröst Bed is for long-term care patients.





#### Domus 4


The Domus 4 Alternating Pressure Redistribution System by Apex is a low-airloss mattress ideal for use in the treatment and prevention of pressure ulcers.

It's advanced, digital multi-mode pump offers alternating pressure redistribution and a static mode with constant low pressure reduction. The Domus 4 offers advanced features such as seat inflation and max firm settings.



### Atlas Hybrid

The Atlas is a non-powered, self-adjusting mattress that utilizes equalization technology to provide maximum pressure redistribution for the prevention and treatment of pressure injuries.

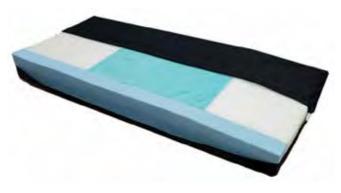


#### Serene Air

Serene Air is a dual mode system with Alternating Pressure Therapy (APT) or Pulsation options. APT modes assist in wound healing and prevention and there are four operating cycles (5/10/15/20 min.) based on therapy requirements. The Pulsation mode increases and decreases the surface's air volume every 15 seconds. Serene Air also has a Constant Low Pressure (static) option, as well as a max inflate and seat inflation for the Fowler position. It also includes a zippered full-length pocket to include a 1" high foam base under the cells.





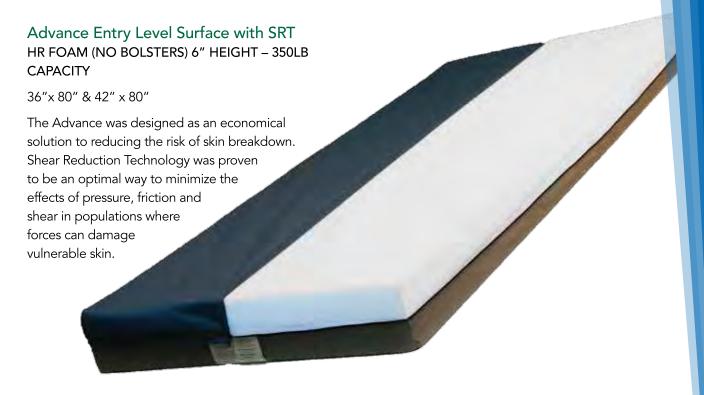

### Adapt Adaptable surface with SRT

3 REMOVABLE GELFLEX SECTIONS INCLUDED - 350LB CAPACITY

36"x 80", 42" x 80", 48"x 80", 54" x 80", 60"x 80



## Bariatric High Density foam surface with SRT 600LB CAPACITY

42" x 80"

Align Bariatric Mattress with Shear Reduction Technology

The Align Bariatric reactive foam replacement surface with Shear Reduction Technology (SRT) is an innovative solution designed to address shearing forces present during client movement in bed or when changing the position of the bed frame, at the same time ensuring proper pressure redistribution. The Align Bariatric is an excellent choice for the heavier client at risk of experiencing skin breakdown. The side bolsters complete with accordion cuts provide an excellent transfer edge. The surface is designed to absorb forces that damage skin, thereby reducing the risk for skin breakdown.

## GelFlex Overlay 2"x36"x80" GelFlex topper covered with Recovery 5 350LB CAPACITY

Our 2"x36"x80" overlay is a perfect alternative to replacing an entire surface or as an added layer of protection for that at-risk individual. Completely covered with Recovery 5 it's a maintenance free option that provides excellent pressure management, comfort and the cooling benefit of gel. Shipped rolled in a box it's easy to store on a shelf and bring out at a moments notice.

### Fall Mat High Density Foam

COVERED WITH VINYL, NON SLIP BOTTOM, FOLDABLE

1"x32"x70"

High-density foam mats are the perfect addition to any bedside. They provide an extra layer of protection in the prevention of injury. The chamfered edge allows for wheelchairs, rollators or lifts to be easily rolled onto the surface.

A rubber base eliminates potential sliding. They are foldable for easy storage.

#### Bed Rail Pad 28"

SIDE RAIL PAD COVERED WITH RECOVERY 5, PADDED WITH GELFLEX

12" x 28"

Our 28" side rail covers are padded on the inside and top edge with GelFlex visco-elastic foam and covered with Recovery 5. Easy to clean and designed to fit most manufactures of beds these pads provide that extra layer of protection in the prevention of injury. Am 8" strap is provided to securely attach them to the bed rail.



#### **Permobil Surfaces**

#### **CUSTOM SIZES**

The standard mattress size is  $36" \times 80"$  for all 5MM products. There is a wider version of  $39" \times 80"$  for both the Adapt 2 and Original Snap Series.

However, for all 5MM matresses, custom sizes are available. Please call for a quote.

## RECOVERY 5 INCONTINENT COVERS AND COVERLETS

All of our mattresses come standard with a Recovery 5 anti-microbial base cover which provides incontinent protection for up to four hours, and aids in keeping skin cool and dry. It is removable, machine-washable, fluid-proof and fire-retardant. This material helps wick away moisture and increases air flow and breathability.

You can add a coverlet to your mattress which has stretchy corners and fits over the mattress (and inserts) like a regular fitted sheet. Each insert can also have its own removable zippered cover.

UPGRADE to the R5 Ultra! The new Recovery 5 Ultra is an additional zippered enclosure cover for the full mattress surface (and inserts) with 4-way stretch, is machine washable, and provides up to 6X more incontinent protection. To order, simply add the code -R5U to the mattress part numbers.

#### **ROHO® MATTRESS SECTIONS**

We have expanded our line of ROHO® compatible mattresses to give you even more options.

These non-powered, adjustable, advanced pressure relieving mattress sections can be inserted into our: Acute, Adapt 2, and Original Snap Series.

These mattress sections are flame and fluid resistant and can be used in a number of configurations – whatever best suits your needs!

#### SPECIFICATIONS AND WARRANTY

Mattress Weight Capacity: 350lb

Standard Mattress Size: 36" x 80",

Custom Sizes Available

Limited Warranty: Two year complete coverage on

mattress and covers

#### TO PLACE YOUR ORDER

Permobil Canada is the exclusive distributor for 5 Minds Mobility's mattresses and

pillows. To order, please contact us at sales.ca@permobil.com.

To learn more, visit our websites at www.permobil.ca and www.5mindsmobility.com.





### Adapt 2 Series Mattress System – 36" x 80"

Size:  $36" \times 80"$ ; Custom sizes available – call for a quote

30" wide insert options: Foam, GelFlex™ or ROHO® (NEW! Foam and GelFlex™ inserts are flippable for longevity and ease of use). A variety of combinations are available to meet your needs using these insert options - just make sure the size adds up to 76"! Note: There is a 4" GelFlex foot bolster at the bottom which adds to the 80" length. It holds the inserts in place without the need for snaps, nor the need to purchase an additional section to complete mattress.

Weight Capacity: 350lb

Cover options: Standard R5 cover provides incontinent protection for up to four hours.

UPGRADE to the R5 Ultra! The new Recovery 5 Ultra top cover for the full mattress surface with 4-way stretch, is machine washable, and provides up to 6X more incontinent protection.



Adapt 2 surface and inserts include a Fire Barrier sock fabric to meet full-scale fire retardant standards (not pictured).

36" base size








Adapt 2 base	Adapt 2 system 1	Adapt 2 system 2	Adapt 2 system 3
AD2-BASE3680	1 - 30" insert	1 - 30" insert	1 - 19" insert
	1 - 8" insert	1 - 8" insert	fits 3 ROHO® inserts
	1 - 19" insert	fits 2 ROHO® inserts	fits 1 ROHO® insert





Adapt 2 system 4	Adapt 2 system 5
fits 4 ROHO® inserts	1 - 30" insert
	1 - 8" insert
	fits 2 ROHO® inserts

Part Number	Description
AD2-BASE3680	36" x 80" x 6" BASE
GELFLEX INSERTS	30" Wide
AD2-830GFS	8" x 30" GELFLEX SECTION
AD2-1930GFS	19" x 30" GELFLEX SECTION
AD2-3030GFS	30" x 30" GELFLEX SECTION
FOAM INSERTS	30" Wide
AD2-830FS	8" x 30" FOAM SECTION
AD2-1930FS	19" x 30" FOAM SECTION
AD2-3030FS	30" x 30" FOAM SECTION
ROHO INSERTS	30" Wide
MATTSECT1610	19" X 30" ROHO® NEOPRENE MATTRESS SECTION
COVERS	
AD2-TC3680	ZIPPERED RECOVERY 5 TOP COVER 36" x 80"
AD2-TC3680-R5U	ZIPPERED RECOVERY 5 ULTRA TOP COVER 36" x 80"
AD2-BC3680	REPLACEMENT RECOVERY 5 BASE COVER 36" x 80"
COV-MATTSECT1610	19" X 30" R5 COVER FOR ROHO® SECTION

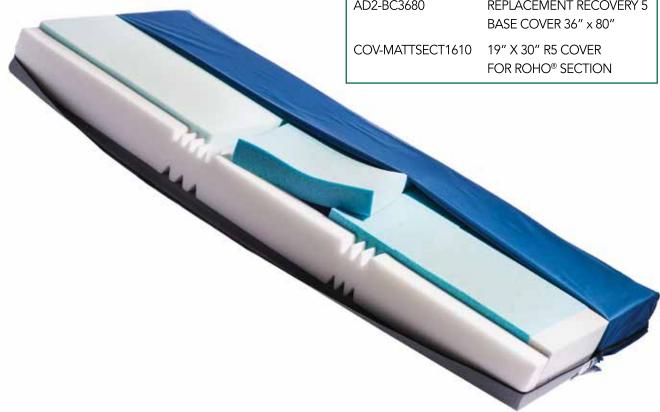

### Adapt 2 Series Mattress System – 39" x 80"

Size:  $39" \times 80"$ ; Custom sizes available – call for a quote

30" wide insert options: Foam, GelFlex™ or ROHO® (NEW! Foam and GelFlex™ inserts are flippable for longevity and ease of use). A variety of combinations are available to meet your needs using these insert options – just make sure the size adds up to 76"!

Note: There is a 4" GelFlex foot bolster at the bottom which adds to the 80" length. It holds the inserts in place without the need for snaps, nor the need to purchase an additional section to complete mattress.

Weight Capacity: 350lb

Cover options: Standard R5 cover provides incontinent protection for up to four hours.

**UPGRADE** to the R5 Ultra! The new Recovery 5 Ultra top cover for the full mattress surface with 4-way stretch, is machine washable, and provides up to 6X more incontinent protection.









Adapt 2 base	Adapt 2 system 1	Adapt 2 system 2	Adapt 2 system 3
AD2-BASE3980	1 - 30" insert	1 - 30" insert	1 - 19" insert
	1 - 8" insert	1 - 8" insert	fits 3 ROHO® inserts
	1 - 19" insert	fits 2 ROHO® inserts	
	fits 1 ROHO® insert		

### Koppla – a Signature Mattress by Permobil

Passive - foam mattress.

Three layers of foam core, 100% medical grade.

- Top layer has 2 inches with a new foam technology (Serene).
- Middle layer has 2 inches.
- Base layer has 2 inches.

Lateral firm boarders.

Top cover: Startex!

Weight capacity: 350lbs.

36 x 90 x 6 inches.

• Other sizes available.



### Sömn – a Signature Mattress by Permobil

Passive - foam mattress.

Two layers of foam core, 100% medical grade.

- Top layer has 2 inches with lower IFD at the heel.
- Base layer has 4 inches.

Lateral firm boarders.

Top cover: Startex!

Weight capacity: 350lbs.

36 x 80 x 6 inches.

• Other sizes available.



### Balanced Aire Non-Powered Self Adjusting Convertible Mattress

- The Balanced Aire mattress replacement is available as either a powered or non-powered support surface.
- The (12) 5.5" tri-laminated foam filled air cells are separated into three distinct pressure zones, each zone independently self-adjusting to the patient's needs.
- The air cells are contained within a high density foam perimeter which creates a smooth transition from the side rail to the air cell.
- In addition, the system is sloped downward in the foot section to reduce pressure in the heel zone.
- Add an alternating pressure pump and the Balanced Aire easily converts to a dynamic powered system for additional pressure redistribution.
- Meets Federal Fire Code CFR 16 part 1633.



This Non-powered Self-adjusting E0373 Group II Surface is Easily Convertible to a Powered E0277 APM

Item#	Description	Dimensions	HCPCS	UOM
BA9600-NP	Non-Powered Self Adj. Air/Foam Mattress	80" (L) × 35" (W) × 7" (W)	E0373	1/ea
BA9600-NP-42 42"	Non-Powered Self Adj. Air/Foam Mattress	80" (L) × 42" (W) × 7" (W)	E0373	1/ea
BA9600-NP-48 48"	Non-Powered Self Adj. Air/Foam Mattress	80" (L) × 48" (W) × 7" (W)	E0373	1/ea
BA9600-NP-54 54"	Non-Powered Self Adj. Air/Foam Mattress	80" (L) x 54" (W) x 7" (W)	E0373	1/ea
BA9600-P	Powered Alternating Pressure Air/Foam Mattress	80" (L) × 35" (W) × 7" (W)	E0277	1/ea
BA9600-P-42 42"	Powered Alternating Pressure Air/Foam Mattress	80" (L) × 42" (W) × 7" (W)	E0277	1/ea
BA9600-P-48 48"	Powered Alternating Pressure Air/Foam Mattress	80" (L) × 48" (W) × 7" (W)	E0277	1/ea
BA9600-P-54 54"	Powered Alternating Pressure Air/Foam Mattress	80" (L) x 54" (W) x 7" (W)	E0277	1/ea

Specifications		
Pump Power	When Used w/Power Source (P)	110 VAC, 60 Hz
Pump Airflow	When Used w/Power Source (P)	5 LPM
Pump Cycle Time	When Used w/Power Source (P) 10 N	<i>f</i> linutes
Pump Type	When Used w/Power Source (P) Compressor	
Pump Alarms	N/A	
Mattress Cover Material	Fluid-Resistant Stretch	
Weight Capacity	BA9600-NP 500lbs.   BA9600-NP-42 600lbs.   BA9600-NP-48 600lbs.   BA9600-NP-54 650lbs.	
Weight Capacity	BA9600-P 500lbs.   BA9600-P-42 600lbs.   BA9600-P-48 600lbs.   BA9600-P-54 650lbs.	
Warranty	5 Year Limited Mattress   1 Year Pum	р



#### **Privacy Screen**

Manufactured with flame retardant bacteriostatic white vinyl and sturdy 1" anodized aluminum tubing, the 3-panel privacy screen from Drive Medical is well-built and functional. The materials are as easy to maintain as they are useful. Each hinged panel measures 56" high and 27.5" wide, giving the product an overall height of 70.5". Four 3" hooded casters make transportation and repositioning a cinch.

#### **Features**

- Each hinged panel measures: 56" high 27.5" wide
- Easy to maintain and clean
- Four, 3" hooded casters
- Overall Height: 70.5"
- Screens are 6 mil flame retardant bacteriostatic white vinyl
- Sturdy, 1" diameter anodized aluminum tubing



#### **Drive Overbed Table**

The large, walnut, wood-grain low pressure laminate top provides a roomy surface for this non-tilt overbed table by Drive Medical. The base of the product is highly adjustable; it can be raised with just the slightest upward pressure and adjusts anywhere from 28" to 45".

The 15" x 30" table top locks in place securely the moment the height adjustment handle is released. The strong, steel "H" base and mast are available in chrome or silver vein finishes. On both models, 2" swivel casters allow for easy maneuverability.



### **POSITIONING CHAIRS**

#### Comfort is Essential For Quality of Life.

BRODA's unique line of tilt-in-space positioning chairs have been improving people's lives for more than 35 years. We have witnessed the positive impact BRODA chairs have had in countless patients' lives. For some, having the appropriate BRODA chair means relief from chronic pain, for others, it's maintaining independence, and in many cases, it means getting out of bed for the first time in months to engage in their surroundings.

Regardless of age, patients seated for extended periods of time deserve to be comfortable. Comfort reduces stress, improves recovery, well-being, and is therefore... absolutely essential to one's quality of life.

Comfort is the key, and when it comes to providing comfort, BRODA is essential.



#### Revive Bariatric Shower Commode

Durable stainless-steel frame and components makes the Revive Bariatric ideal for use in shower rooms. With 14° of posterior tilt, 8° of anterior tilt, up to 65° of recline, full elevating legrest, and removable rotational armrests, the Revive Bariatric provides positioning for safe, thorough bathing, improved ability to wash hair and improved access to abdominal and peri areas.

The padded commode seat and height adjustable headrest provide enhanced patient comfort.

The Revive Bariatric provides an 800lb weight





### Tranquille Auto-Locking Glider

The Tranquille Auto-Locking Glider provides safety for individuals requiring long-term care. The Tranquille provides a gentle gliding action to help create a calming, soothing effect for agitated individuals. The gliding action is excellent for Alzheimer's and dementia patients as it helps reduce the need for wandering.

Designed for fall prevention, the Tranquille Auto-Locking Glider begins to glide once the resident sits in the chair, and stops gliding once the resident begins to vacate the chair. It also features a caregiver lock to prevent unwanted gliding.



### **PATIENT ROOM**

### **RAZAP**

### (Also available with 600lb weight capacity RAZ\_AP600)

The Raz-AP mobile shower commode chair comes standard with 4" of height adjustment, a Molded Harmony Back, 5" dual-locking casters, flip-up padded armrests, and a commode pan. The easy-to-clean, polished stainless steel frame is designed for easy side and front access by the user or attendant.

• Stainless steel frame and hardware



### **PATIENT ROOM**

#### **RAZ AT-Tilt**

#### (Also available in 600lb weight capacity RAZ\_AT600)

The Raz-AT (Attendant Tilt) mobile shower commode chair tilts up to 40° to provide extra pressure reduction and comfort. The unique frame design utilizes locking gas springs that help lift the seat to an upright position, making the Raz-AT one of the easiest tilting mobile shower chairs to operate.

- 18" base and seat frame width of so that the chair can still be positioned over a toilet.
- Comes standard with Raz's exclusive Ischial & Pelvic Alignment System (IPAS), provides 2" (5cm) of fore/ aft 1" (2.5cm) of lateral seat adjustment
- Easy-to-clean, polished stainless steel frame fits over most toilets and is designed for easy side and front access by the user or attendant.
- Harmony contoured molded back is standard
- 4" of height adjustment
- Locking gas struts to facilitate tilting
- 5" Dual-locking casters



#### **EZPZ-AT**

The EZPZ- AT (Attendant Tilt) mobile shower commode chair offers exceptional comfort and function at an affordable price. The EZPZ-AT mobile shower commode chair tilts up to 35° to provide pressure redistribution, a contoured composite seat, and a corrosion-resistant aluminum frame. There are many quality options and accessories available to meet individual needs.



#### EZPZ-SP

The EZPZ- AP (Attendant Propel) and the EZPZ- SP (Self Propel) mobile shower commode chairs offer exceptional comfort and function at an affordable price. The easy-to-clean, corrosion-resistant aluminum frame is designed for easy side and front access by the user or attendant.



## **PATIENT ROOM**

### Bariatric Aluminum Rehab Shower Commode Chair

#### WITH FOUR LOCKING CASTERS

- Can be positioned over a standard toilet or used as a portable self-contained commode.
- Floor to seat height can be adjusted 4" without tools.
- White aluminum frame is attractive and rust resistant.
- Padded seat back and arms.
- Height adjustable.

• Comes with swing-away footrests.	Overall Width
	Overall Height
	Seat to Floor Height
	Seat Dimensions
	Backrest Height
4	Net Weight
	Gross Weight
	Weight Capacity
	Warranty



Accessories		
RTL12085	Biodegradable Sanitary	12 liners/box,
	Commode Liner	6bx/cs



#### Molift Quick Raiser 2

Sit-to-stand hoist with electrically adjustable leg base.

#### Stable and easy to manoeuvre.

Molift Quick Raiser 2 with electrically adjustable leg base is compact, stable and easy to manoeuvre. The unique wheel position and near perfect weight distribution of the base provides a good turning radius and easy manoeuvrability even in narrow spaces.

#### Fits almost everywhere.

The inclined straight column imitates the natural pattern of movement and fullfils the hoisting needs of users with some degree of weight bearing capacity. Molift Quick Raiser 2 has a safe working load of 160 kg (352lbs).

The low base, 10 cm (4"), fits under most beds and chairs while the electrically adjustable legs enables accessibility almost everywhere.

### Hoyer Advanced Power Floor Lift

The Advance™ patient lift from Hoyer represents a true step forward in patient lift design and performance. It uses the very latest design and manufacturing technologies to ensure the Advance is the most functional and stylish lift in its class today.





#### Savaria FL Fixed Lift

A lithium-ion battery powers superior performance

- More lift cycles: delivers 50% more lift cycles per charge than industry average
- Less charging time: recharges from full depletion in only 2 hours, with daily charge needed of only 15 to 30 minutes
- Fewer replacements: battery replacement estimated at every 5 years with normal usage
- More reliable: fewer service calls for dead batteries
- Greener: no lead acid disposal issues and longer life of battery
- No memory effect: battery does not deteriorate from constant charging
- Power: fast lift speed even at maximum capacity lifting weight

With fewer service calls and far fewer battery replacements, the Savaria ceiling lift easily offers lower cost.

#### No-fray lift strap

Savaria ceiling lifts feature a no-fray strap design with channel stitching, preventing accidental twisting of the strap and reducing abrasion found in traditional designs. The long lasting lift strap eliminates one of the most common early life service issues of replacing a frayed strap.

#### Waterproof handset

Easy to clean and waterproof (IP67), the handset is designed with hospital standards in mind. A robust membrane cover prevents moisture ingress and the durable design is built for long life.

#### Easy to service

From installation to required inspections or a service need, the ultra-lightweight design makes everything easier and faster. Servicing can be done on a fixed lift without removal from the track, often without any tools.

The Savaria portable ceiling lift was designed from

the ground up to be a truly portable lift - one that can be easily carried from room to room or place to place when needed. At the heart of the system, a remarkably lightweight lithium battery system delivering ultimate patient lifting power.

### Hoyer® Calibre®

The Hoyer® Calibre® is a uniquely designed electronically operated patient lift for the safe lifting of heavier patients. It is both sturdy and robust with a Safe Working Load rating of 850 lbs.

#### **Features**

- Wrap-around handle allows for multiple caregiver access to maneuver
- Smart Monitor technology drives user compliance and lift longevity
- Intuitive Push Pad assists with initiating movement
- Swan Neck style legs allow Calibre to get closer to the resident
- Weight scale comes standard on Calibre



### Hoyer® Presence®

The Hoyer® Presence® lift has been engineered and designed to handle almost any resident-handling task. Included is Smart Monitor technology that displays lifting activity as well as reminds you to recommended service.

The 500 lb. Safe Working Load and greater spatial area make this lift a perfect fit for larger residents. The Presence comes standard with a powered base and is available with a scale option.

#### **Features**

- Lightweight aluminum construction
- Smart Monitor technology drives user compliance and lift longevity
- Intuitive Push Pad assists with initiating movement.
- Swan Neck style legs allow Presence to get closer to the resident
- Low Leg Option Kit for the Hoyer Presence lift is designed to fit under beds, stretchers, and other furnishings with low clearance



### Hoyer® Stature®

The Hoyer® Stature® combines the benefits of Vertical Lifting and "Clip" style slings. Combined with the ergonomic and intuitive features of the Hoyer Pro family of floor lifts, the Stature is reagrded as the best value in it's class. With a 500 lb. Safe Working Load and one of the largest lifting ranges available, the Stature is able to cope with even the most demanding and technical of resident handling situations. The Stature can lift residents from the floor as well as accommodate transfers to higher surfaces, while allowing the choice of different sling systems.



### One Way Glide

One Way Glide, tube non-slip 45x45cm 17.7x17.7in. EIM38/8

One Way Glide Long 30x130cm 11.8x51.2in EIM30/8

One Way Glide Long 40x130cm 15.7x51.2in. EIM34/8

### Single Use Guide

4Way Glide LPL Nylon Sheet w/NS-lock NS-sides 85x200cm 33.5x78.7in. EIM85/200LPLNS



### **Support Belts**

Support Belt XS, 4 handles, waist size 65-105cm Support Belt S, 5 handles, waist size 65-105cm Support Belt M, 6 handles, waist size 90-130 Support Belt L, 8 handles, waist size 125-150cm Support Belt XL, 11 handles, waist size 145-170cm



#### **Etac Turner PRO**

A sit to stand aid with functional design.

Etac Turner PRO is a turn aid with a functional design that offers safe patient turning with standing support. This turn aid enables a good posture and minimises the risk of injury for the carer when transferring, moreover, it encourages natural standing and sitting movement for the user. The only adjustment necessary

is the height of the leg support and the handle. Thereafter, the transfer can be performed safely whilst the caregiver and user can maintain eye contact and communication.

The orange handle provides a strong, high contrasting colour which can help users with impaired vision or dementia.

The small footplate makes it the perfect choice of transfer aid when transferring in tight spaces. Etac Turner PRO is easy to transport due to its low weight of only 7,5 kg, furthermore it is

compact and easy to grip.

#### Molift Raiser

Transfer platform for active transfer. Safe and easy to use.

Molift Raiser platform is designed to safely move the user from a sitting to standing position while maintaining the best possible ergonomic moving and handling conditions for the caregiver. The unique leg support is designed to be easily adjusted with one hand. Similarly, wheel locks (right or left) will lock both wheels simultaneously providing a safe braking solution in space limited environments. The fine grain finish of the footplate provides a non-slip surface and the handle offers many different grip possibilities. The Molift Raiser is suitable for use with a single carer and our User Guide has been updated to reflect this. These features enable the user to feel safe throughout the entire transfer.



### **BATHING**

#### **Shower Panels**

The CPC1000 console has a thermostatic mixing valve for safe operation without risk of burning. Optionally, it is possible to add a spray disinfection including a 4L reservoir with lock and a collector (sink drain) container for proper use with a shower stretcher.



The CPC2000 console has a thermostatic mixing valve for safe operation without the risk of burning, a display of the temperature of the water and a shower integrated disinfection including a disinfectant tank with lock. Optionally, it is possible

to add a porcelain toilet with flush, freeing up space in the bathroom. It is also possible to use with a commode or to perform drainage stretcher shower.



### **Shower Trolley**

The CDA4000 is the most widespread and most standard model stretcher shower, available in electric version for height adjustment and tilt. The stretcher shower has a comfortable mattress, two inches thick and semi-rigid sides remain in place. The maximum capacity is 400 lb and the tray is available in widths of 24 or 28 inches.



### Genius – Shower Trolley

The Genius is a height adjustable shower trolley that helps to provide hygiene care for totally dependent residents. The Genius can be provided in either manual hydraulic or a battery powered version. The Genius accommodates side transfers from a bed or wheelchair to the shower trolley and back again. Included as standard are a soft mattress, side supports, and leveling and tilting device. A shower panels is also available as an option. (Manual version shown.)





#### CM2000

Made in Canada by ConfortMedic, the CM2000 bathtub hoist is made for the side access to the tub. The CM2000 is made of stainless steel and aluminum for easy disinfection and maximum durability. The seat and backrest are molded in a closed-cell resin, resistant, easy to disinfect and soft to the skin. They can be easily removed for cleaning or replacement.



Made in Canada by , the CM1000 bathtub hoist is made for the front access to the tub. The CM1000 is made of stainless steel and aluminum for easy disinfection and maximum durability. The seat and backrest are molded in a closed-cell resin, resistant, easy to disinfect and soft to the skin. They can be easily removed for cleaning or replacement.





### **Authentic**



### Oslo Collection

LONG TERM CARE FURNITURE | BEDSIDE CABINETS AND CHESTS

#### **Bedside Cabinets**

Dimensions: 21"(W) x 18"(D) x 29.5"(H)

ltem #	Description
BDL100A030-A	1-Door, Cherry
BDL100B030-A	1-Door, Oklahoma Oak
BDL300A010-A	3-Drawer, Cherry
BDL300B010-A	3-Drawer, Oaklahoma Oak



Dimensions: 3-Drawer 29"(W)  $\times$  18"(D)  $\times$  29.5"(H) 4-Drawer 29"(W)  $\times$  18" (D)  $\times$  37"(H)

ltem #	Description
CTL300A010-A	3 Drawer Chest, Cherry
CTL300B010-A	3 Drawer Chest, Oklahoma Oak
CTL400A010-A	4-Drawer Chest, Cherry
CTL400B010-A	4-Drawer Chest, Oaklahoma Oak

## Available laminate finishes

Oklahoma Oak













### Oslo Collection LONG TERM CARE FURNITURE | WARDROBES

#### Dimensions:

2-Door 30"(W) x 23"(D) x 78"(H)

1-Door, 2-Drawer 22"(W) x 23"(D) x 78"(H)

ltem #	Description
WDL2L0A020-A	2-Door, Cherry
WDL2L0B020-A	2-Door, Oklahoma Oak
WDL1S2A030-A	1-Door, 2-Drawer, Cherry
WDL1S2B030-A	1-Door, 2-Drawer, Oklahoma Oak



#### Dimensions:

2-Door, 2-Drawer 30"(W) x 23"(D) x 78"(H)

2-Door, 2-Drawer 38"(W) x 25"(D) x 78"(H

Oversized

Item #	Description
WDL2L2A030-A	2-Door, 2 Drawer Cherry
WDL2L2B030-A	2-Door, 2 Drawer Oklahoma Oak
WDL2O2A030-A	2-Door, 2 Drawer Cherry (oversized)
WDL2O2B030-A	2-Door, 2 Drawer Oklahoma Oak (oversized







### Oslo Collection

### LONG TERM CARE FURNITURE | REPLACEMENT HANDLES AND UNIVERSAL HEAD AND FOOTBOARDS

### Replacement Handles

ltem #	Description
PA-AS C010	Handle for Oslo Cabinet
PT-AS W100	Handle for Oslo Wardrobe

### Universal Head and Foot Boards (requires drilling)

ltem#	Description
HFL360A000	H/F Boards, 36", Cherry
HFL360B000	H/F Boards, 36", Oklahoma Oak
HFL390A000	H/F Boards, 39", Cherry
HFL390B000	H/F Boards, 39", Oklahoma Oak
HFL420A000	H/F Boards, 42", Cherry
HFL420B000	H/F Boards, 42", Oklahoma Oak



# Drive DeVilbiss Healthcare Furniture Suites LONG TERM CARE FURNITURE

Standard features of Brownlee Woods, Smokey Hollow, Mill Creek, Brier Hill, Bronze Lake:

- Choice of 6 most popular wood grained finishes compliment most any décor.
- Drawer slides are epoxy coated and have a 100lb. rating.
- Four sided white vinyl wrapped drawer boxes makes it easier for resident to see and identify personal effects.
- Ventilated backs facilitate air circulation to help prevent bacteria growth and reduce odor.
- Fixed Drawer/Door Handles and Knobs are beautifully crafted and easy to grasp.

Available laminate finishes Kensington Maple



Solar Oak



Wild Cherry



Montana Walnut



Figured Mahogany



Windsor Mahogany





### Drive DeVilbiss Healthcare Furniture Suites LONG TERM CARE FURNITURE

#### STANDARD FEATURES

Ventilated Back



Standard Drawer Slide



Euro Door Hinge



Interior Finish



Arched Toe Kick




Flat Toe Kick



**OPTIONS** Upgraded Drawer Slide



Drawer Lock

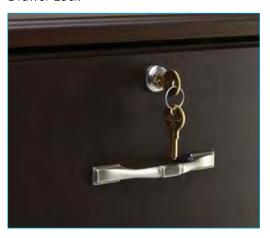

DOOR AND DRAWER HARDWARE

**Knob Styles** 

Satin Nickel | 7476



Antique Brass | 7520



Banded Pewter | 7475



Satin Nickel | 7488



Antique Brass | 7391



Brushed Rope | 7521





#### Royal Ease

How does a caregiver easily move a seated person to-and-from the table? Easily move a seated person to-and-away from the table.

### **Mobility Assistance**

The Royal EZ enables a caregiver to gently move a seated person to and away from the table.

#### Easy to Operate

The Royal EZ can easily move 300lbs with little effort. So easy, a child can do it!

#### **Respect & Dignity**

Since the Royal EZ is very discreet, there is no need for extra attention on guest at meal time.

#### Reduces Injuries

Reduce the risk of injury to staff or caregivers since little effort is required to move a seated person.

#### No Floor damage

Since chairs aren't being dragged, floor and carpet damage is eliminated.

#### Secure & Tough

Many measures were taken during construction and assembly of the Royal EZ to ensure maximum safety.

#### Safety Rear Wheels

The rear wheels are adjusted for tension to prevent chair from moving while not in use.



#### Titan 841SWL-22"-REZ

How does a caregiver assist a person requiring the use of a walker get seated at a table? It Moves. It Swivels. For Individuals with Limited Mobility.

#### It Moves

With the Royal EZ a seated person can be moved to the table with ease.

#### It Swivels

The seat/back of the Titan chair swivels  $90^{\circ}$  left  $/90^{\circ}$  right.

#### 3-Lock Position

The seat locks (centre, 90° left, 90° right), allowing persons to hold onto the arms of the chair as they transfer into the seat

#### Cold Cured Foam

The cold cured foam (CCF) seat cushion prevents water and other liquids from penetrating the cushion.

#### Easy Cleaning

Poly open arms caps allow for easy grip, and a large wipe-out design allows for easy maintenance and clean up.





#### Titan 841CAS-22"

#### TRANSPORT/MOBILITY AID CHAIR

Width	22"
Seat width (between arms)	20"
Seat height	18" or 20"
Arm height	25.5"
Depth	24"
Height	34.5"
Weight	43lbs
Frame rating	800lbs



### TitanPlus+ Equals Space

#### Extra wide seats.

The TITAN Plus+ Series accommodates an even larger guest, with widths of 24–30 inches and a rating of 600lbs. It can be upholstered in any designer fabric and vinyl.

The Titan Plus+ Series embodies the same design characteristics of the Titan Series. With chair widths ranging from 24"- 30", the Titan Plus+ provides the strength, support & comfort required by many of today's full sized guests.

The Titan Plus+ is also available with a complete set of five inch (5") swivel casters, thereby minimizing the effort required on behalf of caregivers to relocate or reposition guests. Each caster can be locked independently by the caregiver.

#### TitanPlus+ Series of Chairs

### Designed for bariatric applications.

The Titan Plus+ Series embodies the same design characteristics of the Titan Series.

With chair widths ranging from 24"- 30", the Titan Plus+ provides the strength, support & comfort required by many of today's

full-sized quests.

Width	33"
Seat width (between arms)	30"
Seat height	18" or 20"
Arm height	25.5"
Depth	24"
Height	35"
Weight	68lbs
Frame rating	600lbs





### Aloe Sleeper Chair

The Aloe Sleeper Chairs effortlessly convert from a chair into a bed. With no mechanical parts, its intuitive design and easy handling means anyone can open and close it without instruction or assistance. Simple to clean, Aloe Sleeper Chairs are constructed from one solid piece of metal, reducing the need for maintenance and repairs.

The Aloe Sleeper Chairs are a superior product that appeal to Infection Control, Environmental Services and Facilities professionals employed in the healthcare continuum.

The Aloe Sleeper Club Chair version combines the ease of use and intelligent design of the Aloe Sleeper Chair with stylish arms that create a relaxing and warm healthcare atmosphere.

A favorite with healthcare staff, the Aloe Sleeper Chairs effortlessly convert from a chair into a bed. With no mechanical parts, its intuitive design and easy handling means anyone can open and close it without instruction or assistance. Simple to clean, Aloe Sleeper Chairs are constructed from one solid piece of metal, reducing the need for maintenance and repairs.

The Aloe Sleeper Chairs are a superior product that appeal to Infection Control, Environmental Services and Facilities professionals employed in the healthcare continuum.



#### **Features**

- Specifically designed for cleanability and infection control.
- A minimalist mechanism comprised of one, stainless steel, highly durable spring.
- Durable one piece frame construction.
- Easy to open and close.
- Locking casters.
- Wall-saver design.
- Replaceable parts.



### **Butterfly Tables**

How do spouses enjoy a meal together when one of them requires the use of wheelchair?



### Chair Caddie

The Chair Caddie enables a spouse/caregiver to easily and safely move a seated person up-to and away-from the table without risk of injury.

- Reduces caregiver fatigue
- Reduces workplace injuries
- Provides dignity to persons with limited mobility

The Chair Caddie enables a caregiver to move a seated person up-to and away-from the table with EASE and SAFETY!



### PromenAid™

PromenAid™ is the first and only modular handrail that can be used anywhere – hallways, stairs, bathrooms, walkways – in dry or wet conditions. Available in a modern champagne-anodized finish, or in a rich architectural bronze finish with antique brass accessories to complement the traditional décor.

PromenAid™ meets the requirements and complies with international residential and commercial standards.









### **■ DISPOSABLES & INFECTION CONTROL**

### Disposable Underpads

This disposable absorbent underpad provides waterproof protection for beds and chairs.



### Quilted Reusable Underpad

- Ideal protection for beds, chairs and other furniture
- Rapidly draws wetness away from the body.
- 100% brushed polyester top layer.
- Extra absorbent polyester and rayon filling, for all levels of incontinence.
- Waterproof backing to prevent leaks.
- Easy to wash and dry.



#### Prevail

First Quality® has become the market leader of Adult Incontinence Products in the North American market by creating value beyond providing products. With a Total Care® Philosophy, First Quality® offers the right products and works with communities to meet both their incontinence program's clinical and financial objectives.

PREVAIL is a full line of innovative product solutions for the lightest to heaviest incontinence needs. Providing high quality product selections that provide comfort, security and dignity for people.

The PREMIER SAMPLE PROGRAM allows one to sample a product prior to purchasing to prevent error in ordering. You can ensure that you are getting the right fit, size and comfort by being able to try the different products available. A guide is available to assist with correct sizing.







#### Evac+Chair

The EVAC+CHAIR® model 300H evacuation chair is a device which permits the safe egress of mobility impaired persons from a building in the event of an emergency or evacuation. The standard model is our lightest weight evacuation chair at only 19lbs. The 300H features a 400lb. weight capacity and a blue powder coated textured finish with a contrasting fire retardant yellow hammock.



### Paraglide

The Paraglide® automated repositioning system for wheelchair users is the latest in assistive device technology that changes the way people are repositioned in wheelchairs in homes, hospitals and care facilities around the world. With Paraglide's patented technology, a wheelchair user can move freely from a slouched position to an upright position with the touch of a button – safely, independently and with dignity.





### DISPOSABLES & INFECTION CONTROL

### Disposable Face Masks, Gloves and Gowns

Harding Medical supplies several types of disposable face masks, gloves and gowns. Surgical masks, gloves and gowns vary by quality and levels of protection. Surgical masks may be labeled as surgical, isolation, dental, or medical procedure masks. Gloves are made of different polymers including latex, nitrile rubber, polyvinyl chloride and neoprene; they come un-powdered, or powdered (to lubricate the gloves), making them easier to put on the hands. Surgical gowns are worn as a barrier to fluid and microbial transmission. They have long sleeves with elastic cuffs and maintain the sterile barrier. Please call for more information to determine what product is best suited for your needs.






### Manorapid

Manorapid is trusted in operating rooms, ICU's, CCU's laboratories and wards in hospitals, clinics and nursing homes across Canada. This key to a great antiseptic not just efficacy but how compatible is it with your skin with repeated use.

Manorapid is used for over 30 years in 70+ countries and plays a key role in controlling the spread bacteria and germs all over the world. Make it your choice to protect your patients and your families.

Manorapid has been an approved product by Health Canada for over 20 years.

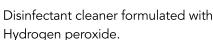

GLDI-H202

Hydrogen peroxide.

Virucidal, fungicidal, bactericidal, kills mold & mildew.

For use on: walls, floors, counters, tables and other hard surfaces, stainless steel.

For use in: schools, daycare centres, office buildings, supermarkets, restaurants, hospitals, long term

care facilities.



### EZ11PlusFully Automatic Autoclave

EZ11Plus, designed to meet the most current sterilization standards ANSI/AAMI ST55, helps practitioners meet today's challenging workloads and regulatory requirements.

#### Valueklave 1730

he Valueklave 1730 will satisfy all of your basic sterilization needs at an economical price. The compact size (7" chamber) makes it a perfect primary unit for small offices or a second autoclave for larger practices. The fast cycle time and heated drying delivers a quick turnaround for wrapped



or unwrapped instruments. The simple controls and easy closing door make this unit a breeze to operate. High speed and reliability are built into every Valueklave 1730.

### Warming Cabinets

Warming Cabinets provide a constant preset temperature within each chamber for surgical

solutions, blankets and other clinical items.



#### Carts

Our durable carts are equipped with a 5-legged base that provides 3 points of constant contact when crossing elevators and thresholds, using long lasting 4" twin-wheel casters. They are ergonomically designed with the end-user in mind, with large amounts of modular space available for accessories.





### Protekt® Finger Pulse Oximeter

#### Features and Benefits

- Ideal for use in multiple applications; hospital, nursing home, community medical treatment, home care, sports center, etc.
- Portable non-invasive spot check oxygen saturation and heart rate.
- Displays digital SpO2 and pulse rate while using bar graph technology to monitor heart rate.
- Automatic power off.
- Two AAA batteries included.
- 1 Year limited warranty.



### 2650KL - Digital Wheelchair Dual Ramp Scale



### 2700KL – Digital Wheelchair Dual Ramp Scale with Extra Large Platform



### **Advantages**

- High Capacity
   High 600lb/270kg capacity provides flexibility to weigh almost any patient.
- EMR Connectivity
   Scale interfaces with a computer, EMR software or other electronic device via various Health o meter® Professional Connectivity Solutions.
- Calculates Body Mass Index (BMI)
   Scale functions include the ability to calculate the patient's BMI.
- Versatile Platform
   Wide and low profile platform accommodates
   wheelchairs and also provides a comfortable
   weighing experience for bariatric patients.
- Value-Added Features
   Operating functions (LB/KG Conversion, lb/KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold/Release, Reweigh, Recall, Auto Zero, Auto Off, Time/Date, Variable Auto Off Time, Audible/Mute Sound Option, Parent/Child Function) included at no additional cost. Digital Wheelchair Dual Ramp Scales 2650KL/2700KL
- Dual Ramps and Remote Display
   Non-slip ramps provide easy on and off access to
   weighing platform and remote display provides
   versatility with the ability to be mounted on a wall
   or placed on a table.
- Everlock® (patent pending)
   Option that allows users to permanently lock the scale in eitherlb or KG to improve safety.



#### Mobilift CX

The Mobilift CX is the easiest way to provide temporary wheelchair accessibility to locations that do not require permanent facilities or where ramps, elevators, and other lifts are not practical/economical. In most cases you can simply roll the Mobilift into position and start providing access immediately – without the headaches of modifying your residence or building.

The Mobilift CX can be used indoors or outdoors and has an unsurpassed 25-year track record of safety and reliability. Over 2500 similar units have been deployed in every kind of environment from the coldest Alaskan winters to the hottest tropical summers. The Mobilift's proven design and durable manufacturing ensure that it will work anytime it is needed.



### HÅG Capisco

The unique shape of the backrest provides perfect support in numerous sitting positions. The saddle seat can be adjusted between low and high working positions, all the way up to standing, enabling you to be more dynamic in your movements. Often referred to as the "Sonographer's Chair," the backrest design

provides full truck and forearm support when seated backwards, and full support for the arm when seated sideways.

#### **Features**

Backrest: Height of backrest (18.1")

Height of foremost support on backrest (5.9" – 9.5")

Seat: Seat depth (13.4" - 16.5"), Seat width (18.5")

HÅG inBalance Mechanism®:

Seat depth

Back height adjustment

Backward tilt tension

Two-position tilt lock

Pneumatic lift:

Standard – 150mm pneumatic lift (16.4" – 21.5")

Options: 200mm pneumatic lift (19.1" - 26")

265mm pneumatic lift (22.4" - 31.9")

Base: Footbase diameter (28.7")

Dual wheel urethane casters

Foot support provides a comfortable way to vary the

position of your feet

Options:

HÅG StepUp® – provides two extra levels for placing

our feet

HÅG Foot Ring – extra level for rest and support for feet Standard fabric color is grade 1 black upholstery or

grade 1 grey vinyl

Color options available in Corner Block fabrics No COM accepted (Customers Own Material)



### Health o meter Portable Digital Chair Scale

#### Advantages

- Easy to Move Lightweight frame with 4 casters make portability easy for bringing the scale to the patient's location.
- High Capacity
   High 600lb/270 kg capacity provides flexibility to weigh almost any patient.
- No Special Shipping Requirements
   The 594KL-c can be shipped via ground service, no need for a pallet or loading dock.
- Portable Digital Chair Scale, Canada Version 594KL-C
- Designed for Patient Comfort
   Large seat with moving armrests allow for easy access on and off the scale.
- Convenient Storage Features
   Folding seat back and armrests make the scale simple to store.
- Value-Added Features
   Operating functions (LB/KG Conversion,lb/KG Lock Out, Body Surface Area (BSA), BMI, Zero, Tare, Hold/Release, Auto Zero, Auto Off) included at no additional cost.

#### **Specifications**

• Capacity: 600lb/270 kg

• Resolution: 0.2lb/0.1 kg

Seat Size (w x d): 18" x 18 ½"
 457mm x 470mm

Product Footprint (w x d x h): 23 ½" x 40 ½" x 37 ½"
 590mm x 1020mm x 960mm

• Product Weight: 39lb/18 kg

• Display: 13/8"/35mm LCD

 Wheels: 4 Casters for ease of movement

 Power Source:
 6 AA batteries (included)

• Functions:lb/KG Conversion, lb/KG Lock Out, BSA, BMI, Zero, Tare, Hold/Release, Auto Zero, Auto Off,

• 2-Year limited warranty



### upCentric Mobile Workstation

The upCentric Mobile Workstation is the perfect solution for people on the move. It allows for flexible movement around a work environment without sacrificing work space. The lockable twin wheel castors keeps the workstation secure when not in motion, and the sealed gas cylinder allows for quick, smooth transition from sitting to standing.

### Features

• Height range: 30"-40"

• Work surface: 28 1/8" x 18 7/8"

• Weight capacity: Up to 88 lbs

• Top and frame colour: Light grey and silver

• Lift: Gas spring for smooth transition s

• Base and casters: One piece cast base, lockable twin wheel casters

Assembly time: 5-10 minutes

• Warranty: 10 years





A division of IMP Group Ltd.



Harding Medical is now offering AIR MILES® Reward Miles.™